

NMR Pipets

Extended tip NMR pipets for retrieving and adding samples in NMR tubes.

Article Number	Wilmad Equivalent	Tip Length	Pkg. Qty.
897085-0007	803	7"	100
897085-0009	803A	9"	100

NMR Sample Tube Rack, 3, 5 and 10 mm

Polypropylene rack with carrying handles features an upper deck with holes and a lower deck support for tube ends. Dimensions: $8\frac{3}{8}$ "L x $4\frac{1}{2}$ "W x $8\frac{3}{4}$ "H.

Article Number	Wilmad Equivalent	Accommodates	Pkg. Qty.
897080-0005	820-A	72— 3 or 5 mm tubes	1
897080-0010	820-B	18—10 mm tubes	1

NMR Tube Washer

Designed to fit standard filtration assemblies available in most labs. The cap is positioned on the 5 mm NMR tube bottom and inserted into the washer as shown below. Aspirator suction provides the motivating force to direct solvent stream into the 5 mm NMR tube. An air stream aids the drying process after solvent reservoir empties. Supplied with one #8 silicone stopper, but without filter flask.


Article Number

897030-0005 5 mm NMR Tube Washer

Accessories

953760-0000 1000 mL Filter Flask
953763-0000 #8 Silicone Stopper, Pkg./ 5

NMR Tube Washer Typical Application

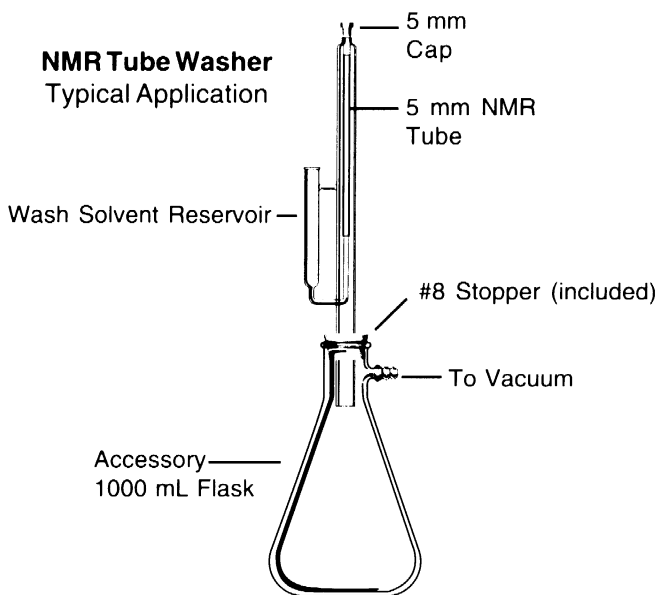

897030-0005